Personenwaage Katharina

Bathroom Scale Katharina | Báscula personal Katharina
Pèse-personne Katharina | Bilancia pesapersona Katharina
Waga osobowa Katharina

Deutsch

English

Español

Français

Italiano

Polski

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BE 1723

Operating instructions | Instrucciones de operación | Mode d'emploi
Istruzioni per l'uso | Instrukcja obsługi

Liebe Kundin, lieber Kunde !

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Dear Customer!

You have chosen to purchase a high-quality product from the ADE brand, which combines intelligent functions with an exceptional design. These bathroom scales make it possible for you to determine and monitor your body weight. The many years of experience of the ADE brand ensure a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

General information

About this manual

This operating manual describes how to safely use and care for your products. Keep this operating manual in a safe place in case you would like to refer to it later. If you pass this item on to someone else, pass on this manual as well.

Failure to observe the recommendations in this operating manual can lead to injuries or damage to your product.

Explanation of symbols

This symbol, combined with the word DANGER, warns against serious injury.

This symbol, combined with the word WARNING, warns against moderate and mild injuries.

This symbol, combined with the word NOTE, warns against material damage.

This symbol signifies additional information and general advice.

Table of Contents

General information	12
Intended Use	14
Safety	14
At a glance	16
Scope of delivery	16
Start-up	17
Weighing	18
Resetting the scale	18
Cleaning	18
Fault/Remedy.....	19
Technical Data	19
Declaration of Conformity.....	20
Warranty	20
Disposal.....	21

Intended Use

The bathroom scale are designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

The maximum carrying capacity of 180 kg must not be exceeded.

Safety

- This item can be used by children 8 years of age and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge, if they are supervised or have been instructed in the safe use of the item and understand the dangers that may arise.
- Children must not play with the item.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries, when swallowed, can be life-threatening. For this reason, the scale and battery should be kept in a place out of reach of children. If a battery is swallowed, medical assistance must be sought immediately.

DANGER for Children

- Keep children away from packaging material. Risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.

RISK of fire/burning and/or explosion

- There is a risk of explosion if batteries are replaced incorrectly. Replace the battery only with one of the same or a similar type. Pay attention to the "Technical Data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Otherwise there is an increased risk of explosion and of leakage; furthermore, gases may escape.

- Avoid contact with skin, eyes, and mucosa. If contact is made with the battery fluid, rinse the affected areas immediately with plenty of clean water and see a doctor straight away.


WARNING against risk of injury

- Do not step onto the outer edge of any one side of the scale. Risk of tipping!
- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!

BEWARE of material damage

- Place the scale on a stable, even surface.
- Do not overload the scale. Pay attention to the “Technical Data”. Place yourself on the scale exclusively for weighing purposes.
- Protect the scale from sharp jolts and vibrations.
- Remove the battery from the scale once it has been used or when you are no longer using the scale. You will thus avoid damage that may result from the leakage of battery acid.
- When inserting and changing batteries, pay attention to the correct polarity (+/-).
- Do not subject the battery to any extreme conditions, by – for example – storing it on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the item. Only allow repairs to be carried out at a specialist workshop. Repairs carried out incorrectly can considerably endanger the user.
- Do not place the scale in the vicinity of devices which may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indicators or malfunctions.
- There are feet on the underside of the scale. Since surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place an anti-slip underlay beneath the scale.

At a glance

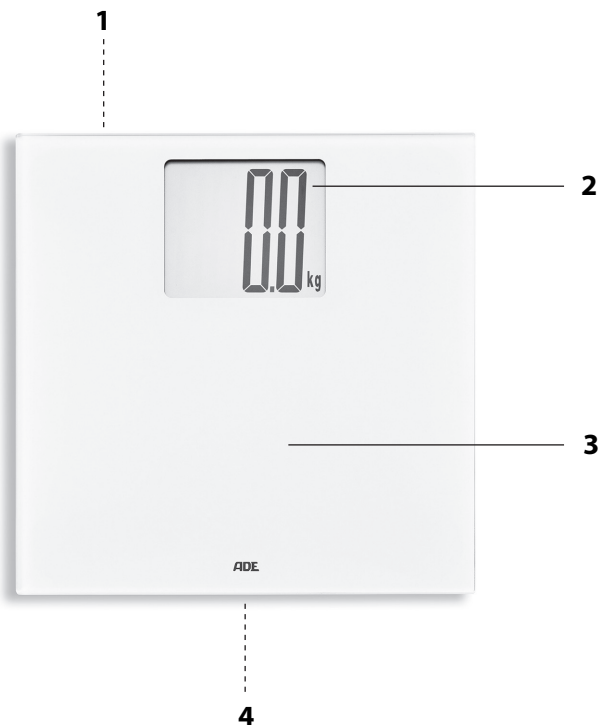- 1 **Unit** button (underside): Changing the unit of measurement
- 2 Display
- 3 Weighing surface
- 4 Battery compartment (underside)

Scope of delivery

- Bathroom scales
- Battery type LR03, 1.5 V  (3x)
- Operating Manual

Start-up

Inserting/changing batteries

For use you require 3 batteries of type LR03 (AAA), 1.5 V (included). The battery compartment is located on the underside of the device.

1. Push the locking tab towards the edge of the scale and remove the battery compartment cover.
2. Insert the batteries as shown on the base of the battery compartment. Pay attention to the correct polarity (+/-).
3. Replace the battery compartment cover. The closing flap must engage audibly.

Suggestions for changing the battery:

- A battery change is necessary when "Lo" or nothing at all is shown on the display.
- Only use the type of battery that is specified in the "Technical Data".
- Always exchange all of the batteries at the same time.
- Dispose of the old battery in an environmentally responsible way (see "Disposal").
- Calibrate the scale, see "Setting the scale".

Choosing the unit of measurement

When the scale is delivered, the unit of measurement is set to "kg" as standard. However, the measurement result can also be displayed in the English unit "lb" (pound) or "st lb" (stone:pound).

1. To switch on the scale, step on the middle of the weighing surface.
2. Press repeatedly the Unit button on the back of the scale.
Note the display.

Positioning the scale

1. Place the scale on an even, firm, and dry surface. Do not place the scale on carpet.
2. First, briefly step onto the middle of the weighing surface to calibrate the scale. As soon as the display shows "0.0", calibration is complete.
This calibration must be performed after each battery change.

Weighing

WARNING against injuries

- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!
- Position yourself on the scale. Remain standing still and do not hold on to anything.

The scale switches on and the display flashes shortly afterwards: The weight measurement is displayed in the selected unit of measure.

A few seconds after you step off the weighing surface, the scale will switch off automatically.

- The weight is displayed for up to approx. 10 seconds after leaving the scale. When you step on the weighing platform using slight force, the display shows "0.0" and you can weigh yourself again.
- Weigh yourself with as little clothing on as possible.
- Try to weigh yourself always at the same time and under consistent conditions in order to obtain comparable values.

Resetting the scale

Reset the scale to factory settings if unusual or obviously incorrect values are shown on the display.

1. Remove the batteries from the battery compartment and wait a moment.
2. Re-insert the batteries.

Pay attention to the correct polarity of the batteries (+/-).

After switching the scale on, you can weigh yourself as usual.

Cleaning

NOTE on material damage

- Do not immerse the scale in water or other liquids.
- Do not under any circumstances use harsh, scraping, or abrasive cleaning agents or the like to clean the scale. If you do, the surface could be scratched.
- Wipe down the scale as needed using a lightly dampened, soft cloth.
- Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

Fault/Remedy

Fault	Possible Cause
No function	<ul style="list-style-type: none"> ▪ Battery empty or inserted the wrong way around?
The scale appear to be displaying false weighing results	<ul style="list-style-type: none"> ▪ Are the scale standing on a firm, even surface? ▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scale? ▪ Remove the battery from the scale and wait a moment. Reinsert the battery.
“Err” is shown on the display	<ul style="list-style-type: none"> ▪ The maximum capacity of the scale (180 kg) has been exceeded.
“Lo” is shown on the display	<ul style="list-style-type: none"> ▪ The batteries are empty (see “Inserting/ changing batteries”).

Technical Data

Model:	Bathroom Scale BE 1723 “Katharina”
Battery:	3x LR03/1.5 V 
Rated current:	2 mA
Measuring range:	max. 180 kg
Variance:	1 % +/- 100g
Indicator division:	100 g
Weight:	approx. 1.6 kg
Dimensions of scale	approx. 300 x 300 x 23 mm
Distributor: (no service address)	Waagen-Schmitt GmbH Hammer Steindamm 27-29 22089 Hamburg, Germany

Our products are constantly being developed and improved.
For this reason, design and technical modifications are possible at any time.

Declaration of Conformity

This scale has been manufactured in accordance with the harmonised European standards.

It complies with the provisions of the EC directives specified below:

- EMC Directive 2014/30/EU
 - RoHS Directive 2011/65/EU
- in the currently valid versions.

This declaration loses its validity if modifications are made to the product without our approval.

Hamburg, January 2018
Waagen-Schmitt GmbH
Hammer Steindamm 27-29
22089 Hamburg, Germany

Warranty

The company Waagen-Schmitt GmbH guarantees for 2 years from the purchase date the cost-free rectification of defects on the basis of material or fabrication errors through repair or exchange.

In the event of a warranty claim, please return the product with the proof of purchase (stipulating the reason for claim) to your dealer.

Disposal

Disposing of the packaging

Dispose of the packaging with similar materials. Put cardboard with waste paper, films with recyclable materials.

Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

Devices must not be disposed of in normal household waste

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Hand in the old device at a collection point for electronic waste or at a recycling depot.

For more information, contact your local waste disposal firm or your local administrative body.

Recyclable and disposable batteries do not belong in household waste

All batteries must be handed in at a collection point in your community, suburb, or at the retailer. Such batteries can thus be disposed of in an environmentally friendly way.

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