## Personenwaage Ylvie

Bathroom Scale Ylvie | Báscula personal Ylvie  
Pèse-personne Ylvie | Bilancia pesapersona Ylvie  
Waga osobowa Ylvie

Deutsch

English

Español

Français

Italiano

Polski

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BE 1722

Operating instructions | Instrucciones de operación | Mode d'emploi  
Istruzioni per l'uso | Instrukcja obsługi

# Liebe Kundin, lieber Kunde !

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# Dear Customer!

You have chosen to purchase a high-quality product from the brand ADE, which combines intelligent functions with an exceptional design. These bathroom scales make it possible for you to determine and monitor your body weight. The many years of experience of the brand ADE ensures a high technical standard and proven quality.

We hope you enjoy weighing yourself and wish you success in reaching your goal weight!

Your ADE Team

## General information

---

### About this manual

This operating manual describes how to safely use and care for your products. Keep this operating manual in a safe place in case you would like to refer to it later. If you pass this item on to someone else, pass on this manual as well.

Failure to observe the recommendations in this operating manual can lead to injuries or damage to your product.

### Explanation of symbols

This symbol, combined with the word DANGER, warns against serious injury.

This symbol, combined with the word WARNING, warns against moderate and mild injuries.

This symbol, combined with the word NOTE, warns against material damage.

This symbol signifies additional information and general advice.

# Table of Contents

---

General information .....	12
Intended Use .....	14
Safety .....	14
Overview .....	16
Scope of Delivery .....	16
Start-up .....	17
Weighing .....	18
Resetting the scale .....	18
Cleaning .....	18
Fault/Remedy .....	19
Technical Data .....	19
Declaration of Conformity .....	20
Guarantee .....	20
Disposal .....	21

## Intended Use

---

The bathroom scale is designed for the weighing of people in the home. It is not suitable for use in clinics, medical practices, etc.

The maximum carrying capacity of 180 kg must not be exceeded.

## Safety

---

- This item can be used by children 8 years of age and above, as well as by persons with reduced physical, sensory, or mental capacities or lacking experience and knowledge, if they are supervised or have been instructed in the safe use of the item and understand the dangers that may arise.
- Children must not play with the item.
- Cleaning and user maintenance must not be carried out by children, unless they are 8 years of age or older and are being supervised.
- Batteries can be life-threatening. Therefore, keep batteries out of reach of small children. If a battery has been swallowed or has entered the body in any other way, seek medical help immediately. Failure to do so within 2 hours may result in severe internal burns that could cause death.

## DANGER for children

- Keep children away from packaging material. Risk of choking if swallowed. Also beware that children do not pull the packaging material over their heads.

## RISK of fire/burning and/or explosion

- There is a risk of explosion if batteries are replaced incorrectly. Replace the battery only with one of the same or a similar type. Pay attention to the "Technical Data".
- Batteries must not be recharged, reactivated by other means, dismantled, thrown into fire, or short-circuited. Otherwise there is an increased risk of explosion and of leakage; furthermore, gases may escape.

- Avoid contact with skin, eyes, and mucosa. If contact is made with the battery fluid, rinse the affected areas immediately with plenty of clean water and see a doctor straight away.

## **WARNING against risk of injury**


- Do not step onto the outer edge of any one side of the scale. Risk of tipping!
- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!

## **BEWARE of material damage**

- Place the scale on a stable, even surface.
- Do not overload the scale. Pay attention to the “Technical Data”. Place yourself on the scale exclusively for weighing purposes.
- Protect the scale from sharp jolts and vibrations.
- Remove the battery from the scale once it has been used or when you are no longer using the scale. You will thus avoid damage that may result from the leakage of battery acid.
- When inserting and changing batteries, pay attention to the correct polarity (+/-).
- Do not subject the battery to any extreme conditions, by – for example – storing it on radiators or under direct sunlight. Increased risk of leakage!
- Clean the battery and device contacts, as needed, before inserting the battery.
- Do not make any modifications to the item. Only allow repairs to be carried out at a specialist workshop. Repairs carried out incorrectly can considerably endanger the user.
- Do not place the scale in the vicinity of devices which may emit electromagnetic radiation, such as mobile phones, wireless devices, or the like. Doing so may lead to incorrect display indicators or malfunctions.
- There are feet on the underside of the scale. Since surfaces are made from various materials and are treated with diverse care products, the possibility cannot be completely excluded that some of these substances contain constituents that attack and soften the feet of the device. If necessary, place an anti-slip underlay beneath the scale.

# Overview

---

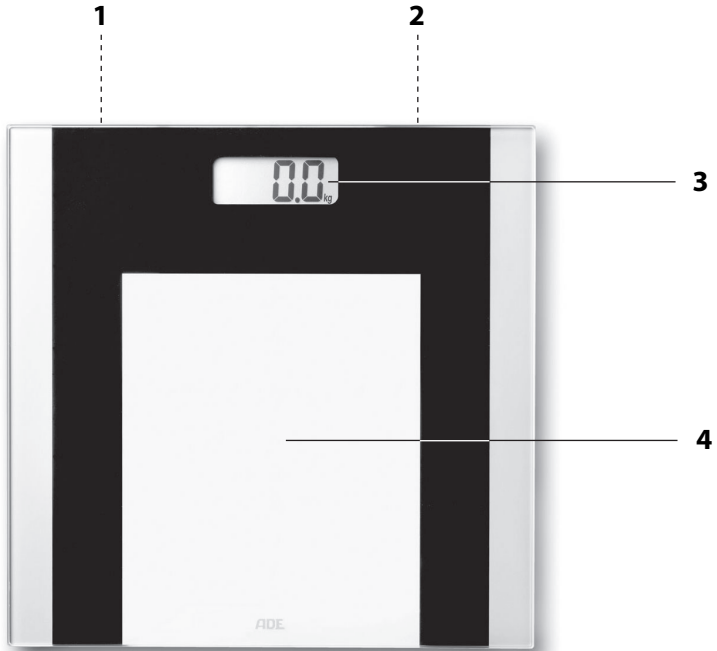- 1 **UNIT** button (underside): Changing the unit of measurement
- 2 Battery compartment (underside)
- 3 Display
- 4 Weighing surface

## Scope of Delivery

---

- Bathroom scales
- Battery type CR2032, 3 V 
- Operating Manual

# Start-up

## Inserting/changing battery

For use you require one battery of the type CR2032, 3 V (included). The battery compartment is located on the underside of the device.

1. Remove the battery compartment cover.
2. Push the battery towards the contact springs under the big, angular retaining plate and then press it down. The battery will snap into the holding lugs on the opposite side and the battery will be held firmly in its holder.

Pay attention to the correct polarity (+/-). The plus pole (+) must point upwards.

3. Put the battery compartment cover back in place. The battery compartment cover must snap securely into place.

Suggestions for changing the battery:

- A battery change is necessary when "Lo" or nothing at all is shown on the display.
- Only use the type of battery that is specified in the "Technical Data".
- Dispose of batteries in an environmentally appropriate way (see "Disposal").
- Calibrate the scale, see "Positioning the scale".

## Choosing the unit of measurement

When the scale is delivered, the unit of measurement is set to "kg" as standard. However, the measurement can also be displayed in the English unit of measure "lb" (pounds).

- Push the **UNIT** button on the back of the scale.  
Pay attention to the display indicator.

## Positioning the scale

1. Place the scale on an even, firm, and dry surface. Do not place the scale on carpet.
2. First, briefly step onto the middle of the weighing surface to calibrate the scale. As soon as the display shows "0.0", calibration is complete.  
This calibration must be performed after each battery change.

# Weighing

---

## **WARNING** against injuries

- Do not step onto the scale with wet feet or when its surface is damp. Risk of slipping!
- Position yourself on the scale. Remain standing still and do not hold on to anything.

The scale switches on and the display flashes shortly afterwards: The weight measurement is displayed in the selected unit of measure.

A few seconds after you step off the weighing surface, the scale will switch off automatically.

- The weight is displayed for up to approx. 10 seconds after leaving the scale. When you step on the weighing platform using slight force, the display shows "0.0" and you can weigh yourself again.
- Weigh yourself with as little clothing on as possible.
- Try to weigh yourself always at the same time and under consistent conditions in order to obtain comparable values.

# Resetting the scale

---

Reset the scale to factory settings if unusual or obviously incorrect values are shown on the display.

1. Remove the battery from the battery compartment and wait a moment.
2. Re-insert the battery. Pay attention to the correct polarity of the batteries (+/-). The plus pole (+) must point upwards.

After switching the scale on, you can weigh yourself as usual.

# Cleaning

---

## **NOTE** on material damage

- Do not immerse the scale in water or other liquids.
- Do not under any circumstances use harsh, scraping, or abrasive cleaning agents or the like to clean the scale. If you do, the surface could be scratched.
- Wipe down the scale as needed using a lightly dampened, soft cloth.
- Disinfect the weighing surface from time to time using a suitable disinfecting agent (available at specialist retailers and pharmacies).

## Fault/Remedy

Fault	Possible Cause
No function	<ul style="list-style-type: none"> <li>▪ Battery empty or inserted the wrong way around?</li> </ul>
The scale appear to be displaying false weighing results	<ul style="list-style-type: none"> <li>▪ Is the scale standing on a firm, even surface?</li> <li>▪ Is a device that can cause electromagnetic disturbances located in the vicinity of the scale?</li> <li>▪ Remove the battery from the scale and wait a moment. Reinsert the battery.</li> </ul>
“Err” is shown on the display	<ul style="list-style-type: none"> <li>▪ The maximum capacity of the scale (180 kg) has been exceeded.</li> </ul>
“Lo” is shown on the display	<ul style="list-style-type: none"> <li>▪ The battery is empty (see “Inserting/ changing the battery”).</li> </ul>
“C” is shown on the display	<ul style="list-style-type: none"> <li>▪ An error occurred. Step off the scale and repeat the weighing procedure.</li> </ul>

## Technical Data

Model:	Bathroom Scale BE 1722 “Ylvie”
Battery:	1x CR2032/3 V 
Rated current:	2.5 mA
Measuring range:	max. 180 kg
Variance:	1 % +/- 100g
Indicator division:	100 g
Weight:	approx. 1.5 kg
Dimensions of scale:	approx. 300 x 300 x 6 mm
Distributor: (no service address)	Waagen-Schmitt GmbH Hammer Steindamm 27-29 22089 Hamburg, Germany

Our products are constantly being developed and improved.  
For this reason, design and technical modifications are possible at any time.

# Declaration of Conformity

---

Waagen-Schmitt GmbH hereby declares that the bathroom scale BE 1722 complies with EU directives 2014/53/EU and 2011/65/EU.

The full text of the EU declaration of conformity is available at the following internet address: [www.ade-germany.de/doc](http://www.ade-germany.de/doc)

The declaration loses its validity if modifications have been made to the device without our approval.

Hamburg, March 2018  
Waagen-Schmitt GmbH

# Guarantee

---

The company Waagen-Schmitt GmbH guarantees for 2 years from the purchase date the cost-free rectification of defects on the basis of material or fabrication errors through repair or exchange.

In the event of a guarantee claim, please return the product with the proof of purchase (stipulating the reason for claim) to your dealer.

# Disposal

---

## Disposing of the packaging

Dispose of the packaging with similar materials. Put cardboard with waste paper, films with recyclable materials.

## Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

### **Devices must not be disposed of in normal household waste**

At the end of its life, the product must be disposed of in an appropriate way. In this way, valuable materials contained in the device will be recycled and the burden on the environment avoided. Hand in the old device at a collection point for electronic waste or at a recycling depot.

For more information, contact your local waste disposal firm or your local administrative body.

### **Recyclable and disposable batteries do not belong in household waste**

All batteries must be handed in at a collection point in your community, suburb, or at the retailer. Such batteries can thus be disposed of in an environmentally friendly way.

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